



4-Hole Rotating Receiver

Instructions For Use

1258



1258



1258-B

1. Application & Description

- This adapter is exclusively for lower extremity prosthetic fittings
- It connects any 4-hole system to a pyramid adapter

2. Technical Specifications

P/N	Material	Build Height	Weight	Weight Limit
1258	Titanium	23.8mm	128g	170kg/375lbs
1258-B	Titanium	23.8mm	128g	170kg/375lbs

- **Recommended Torque:** Always cross tighten alignment set screws
- **M8 Alignment Set Screw Torque:**
 - (Titanium): **15Nm**
- **M6 Flat Head Screw Torque:**
 - (Titanium): **10Nm**
- **M6 Cap Screw Torque:**
 - (Titanium): **10Nm**
- **Thread Locker:** Apply Loctite once alignment and connection are finalized

3. Installation Instructions


- **Angulation:** Adjust angulation alignment accordingly taking caution to loosen adjacent screw to prevent slide galling, which can cause noise after torquing the set screws. Apply Loctite and torque to 15Nm when done.
- **Rotation:** Adjust the desired rotation, tighten the clamp cap screw, and then the 4-hole mounting screws to the recommended torque. Always loosen the two mounting screws on the Clamp Screw side when adjusting rotation, then re-tighten and torque.

4. Safety & Environment

- **Temperature Range:** Between -20°C to 80°C (-4°F to 176°F)
- Rinse with clean water if exposed to salty, chlorinated, or otherwise potentially corrosive liquids
- **Storage:** Store at room temperature; avoid wet; very hot or extremely cold environments.
- **Warranty:** 2 years, provided it is used as intended without modification and only the original user.
- Do not reuse on another prosthesis
- **Conformity:** This product meets the requirements Council EU 2017/745 (MDR) for medical products. This product has been classified as a class I product according to the classification criteria outlined in appendix IX of the guidelines. Please keep this manual in safe place for future reference.



1258 IFU REV.A (5-5-26)
Specifications subject to change.

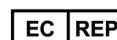

MDSS GmbH