

3030 / 3050 MAGLiners – Transtibial / Transfemoral

HYDRO – Exclusively designed for the MAGLOCK

ST&G revolutionized the liner industry with the use of Parylene® on the surface of the silicone liner. We were also the first to encapsulate the outer cloth into the liner. We have eliminated delamination. We have designed a new cloth specially for this process.

The Hydro MAGLiner's silky smooth mineral oil embedded Platinum cured medical grade silicone liner has the unique ability to repel water. This hydrophobic property makes the Hydro the easiest and fastest liner to clean and sanitize

The Parylene® surface reduces surface friction between the liner/skin and liner/socket for better comfort and donning. The Hydro MAGLiner reduces skin shear for fragile skin and will help with blister reduction. The matrix which is built into the liner prevents distal elongation and delamination.

The Hydro Skin offers proprietary donning and doffing strips to aid the patient in easier application/removal. It is suitable for moderate to high activity patients who require cushioning over sensitive areas in the socket. The Hydro Skin liner's opaque pigment gives a natural "skin tone" to all ethnic skin types.

- 6 Month warranty

Transtibial

P/N	Description	P/N	Description
3030-18	MAGLiner Hydro, 18	3030-28	MAGLiner Hydro, 28
3030-20	MAGLiner Hydro, 20	3030-30	MAGLiner Hydro, 30
3030-22	MAGLiner Hydro, 22	3030-32	MAGLiner Hydro, 32
3030-23.5	MAGLiner Hydro, 23.5	3030-34	MAGLiner Hydro, 34
3030-25	MAGLiner Hydro, 25	3030-36	MAGLiner Hydro, 36
3030-26.5	MAGLiner Hydro, 26.5		

Transfemoral

P/N	Description	P/N	Description
3050 TF - 30	MAGLiner Hydro, 30	3050 TF - 38	MAGLiner Hydro, 38
3050 TF - 32	MAGLiner Hydro, 32	3050 TF - 40	MAGLiner Hydro, 40
3050 TF - 34	MAGLiner Hydro, 34	3050 TF - 45	MAGLiner Hydro, 45
3050 TF - 36	MAGLiner Hydro, 36		

* Using the MAGLINER with the MAGLOCK offers the lowest profile possible.

MAGLiners, PREMIER & HYDRO are designed exclusively for the MAGLOCK Suspension System. The **MAGLOCK** eases the frustration of donning a pin lock liner for patients K1 to K3. The **MAGLOCK** is developed using a rare earth magnet which achieves up to 90lbs of adhesion. The large distal target allows the patient to be 10 degrees off center and still get suspension.



ST&G USA CORPORATION
2691 Saturn Street
Brea, CA 92821

TELEPHONE: (714) 524-0663


FAX: (714) 364-8113

E-MAIL: info@stngco.com

www.stngco.com

3040 / 3060 MAGLiners - Transtibial / Transfemoral

PREMIER – Exclusively designed for the MAGLOCK

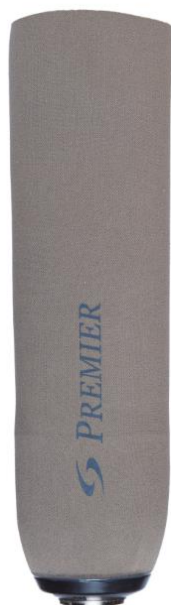

This premium liner is designed for the health of the residual limb and ease of donning. Manufactured with mineral oil embedded Platinum cured medical grade silicone and sheathed in a seamless durable fabric. The matrix which is built into the liner prevents distal elongation.

* 6 Month warranty



Transtibial

P/N	Description	P/N	Description
3040-18	MAGLiner Premier, 18	3040-28	MAGLiner Premier, 28
3040-20	MAGLiner Premier, 20	3040-30	MAGLiner Premier, 30
3040-22	MAGLiner Premier, 22	3040-32	MAGLiner Premier, 32
3040-23.5	MAGLiner Premier, 23.5	3040-34	MAGLiner Premier, 34
3040-25	MAGLiner Premier, 25	3040-36	MAGLiner Premier, 36
3040-26.5	MAGLiner Premier, 26.5		

Transfemoral

P/N	Description	P/N	Description
3060 TF - 30	MAGLiner Hydro, 30	3060 TF - 38	MAGLiner Hydro, 38
3060 TF - 32	MAGLiner Hydro, 32	3060 TF - 40	MAGLiner Hydro, 40
3060 TF - 34	MAGLiner Hydro, 34	3060 TF - 45	MAGLiner Hydro, 45
3060 TF - 36	MAGLiner Hydro, 36		

* Using the **MAGLINER** with the **MAGLOCK** offers the lowest profile possible.

MAGLiners, PREMIER & HYDRO are designed exclusively for the **MAGLOCK Suspension System**. The **MAGLOCK** eases the frustration of donning a pin lock liner for patients K1 to K3. The **MAGLOCK** is developed using a rare earth magnet which achieves up to 90lbs of adhesion. The large distal target allows the patient to be 10 degrees off center and still get suspension.



ST&G USA CORPORATION
2691 Saturn Street
Brea, CA 92821


TELEPHONE: (714) 524-0663

FAX: (714) 364-8113


E-MAIL: info@stngco.com

www.stngco.com


HOW TO MEASURE Transtibial Limb

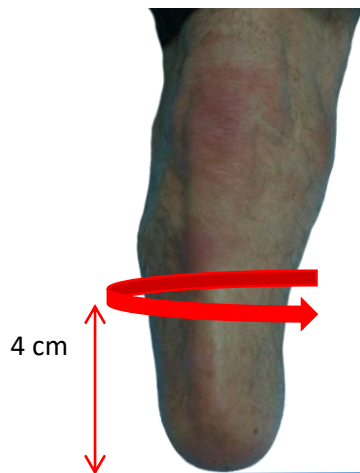


1. Measure the circumference of the residual limb at 4cm from the distal end
2. Choose the liner size, which corresponds to the measurement or one smaller.
3. Never choose a liner size, which is larger than the measurement

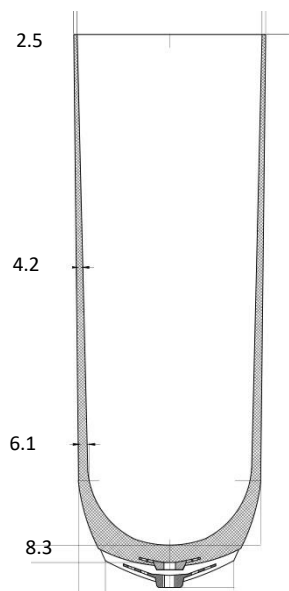


FORTE

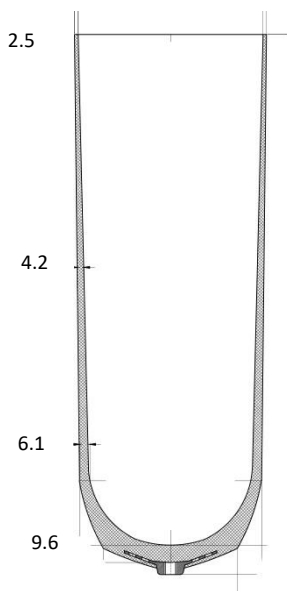

PARASIL

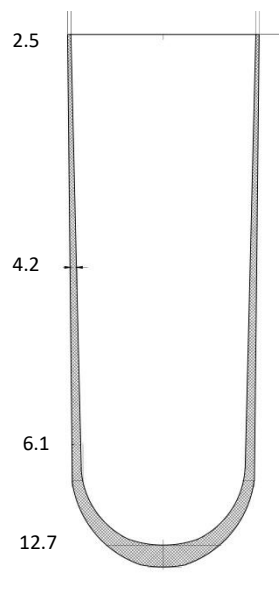

CUSHION

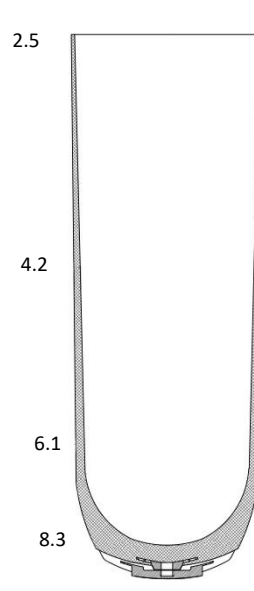

MAGLINER

BK LINERS

WALL AVERAGES


ST&G USA CORPORATION
2691 Saturn Street
Brea, CA 92821

TELEPHONE: (714) 524-0663
FAX: (714) 364-8113
E-MAIL: info@stngco.com




www.stngco.com


HOW TO MEASURE Transfemoral Liner

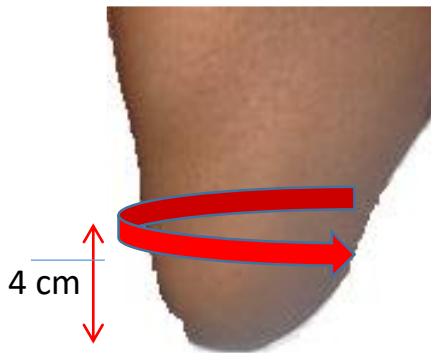


1. Measure the circumference of the residual limb at 4cm from the distal end
2. Choose the liner size, which corresponds to the measurement or one smaller.
3. Never choose a liner size, which is larger than the measurement

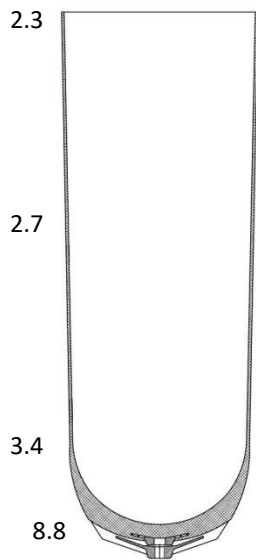


FORTE

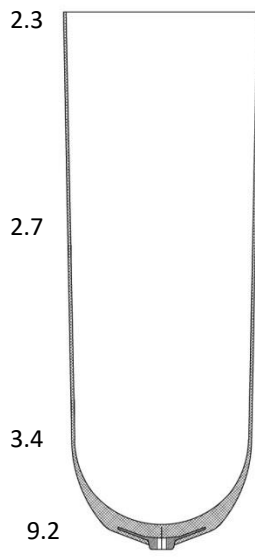

PARASIL

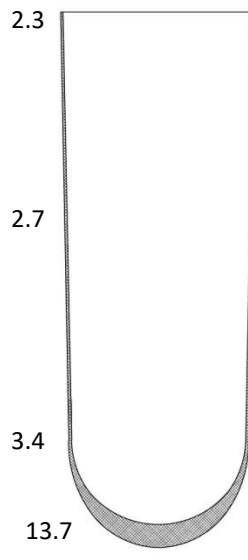

CUSHION

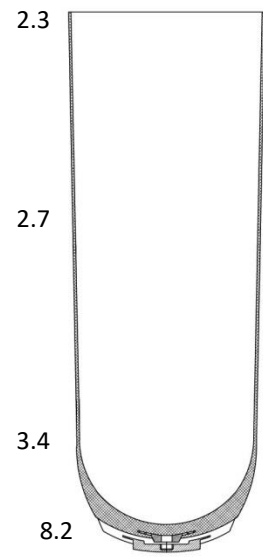

MAGLINER

AK LINERS

WALL AVERAGES

ST&G USA CORPORATION
2691 Saturn Street
Brea, CA 92821

TELEPHONE: (714) 524-0663
FAX: (714) 364-8113
E-MAIL: info@stngco.com

 **ST&G**
CORPORATION
www.stngco.com